

NNMC Bethesda

Medweb

Home

View Patients

Patient schedule

Machine schedule

Monitor

Personal Profile

Statistics

Logout

Administrator Site

Admin Site SSL

Contacts: Dave Cundiff

Business hours: M-F 0800-1600

Announcements: 301.295.5812

Address and driving directions

Address: United States
8901 Wisconsin Ave, Bethesda MD

[Click here to access DICOM studies.](#)

Connect to 131.158.69.236

Patient Record Access

User name:

Password:

OK Cancel

Referring provider
Logs on to
Medweb at NNMC
Bethesda

NNMC Bethesda

[Ophthalmologist]

Medweb tele-medicine intranet software - Version 4.5.20 - Last updated 2005-08-0

New Study

On Call

Preferences

Save settings

View Changes

Manuals

Get latest plugin

Your license has 3380 days remaining.
Your Medweb Server manager can update the license.

Limit: 10	Status: Unread	Find:	Modality: (any)	In last: (any)	Group By: Site/Station	Query: Local
--------------	-------------------	-------	--------------------	-------------------	---------------------------	-----------------

Search: Advanced Current Clear

User: pgunther Drive status: 110057M available for studies
0 images left to send 0 images left to compress

☐ View only pre-selected images

Click on New Study

Name	Date	Time	Mgmt
Site: Usns Comfort Station name: Compact 3 Modality: CR			
<input type="checkbox"/> [Redacted]	CHEST	1	2005-09-10 12:11:06
<input type="checkbox"/> [Redacted]	CHEST	2	2005-09-12 16:38:32
<input type="checkbox"/> [Redacted]	CHEST	1	2005-09-13 11:48:52
Site: Usns Comfort Station name: Compact 2 Modality: CR			
<input type="checkbox"/> [Redacted]	CHEST	1	2005-09-11 12:35:33
Site: Nunc Station name: GEMSOW Modality: MR			
<input type="checkbox"/> [Redacted]	mri, brain w w/o con	139	2005-08-22 08:14:51
Site: Rota Station name: Rota Medweb Modality: SC			
<input type="checkbox"/> Rota Optometry Testimage	1127227696	1	2005-09-20 16:49:44

Click here to access DICOM studies.

ographics Dermatology PATH Ophthalmology Images

ient Billing Military Hx PE

Patient Data

Patient ID 1127401615 Last Name Gunther First Name Peter MI J

DOB 2/3/70 Sex M Marital Status Married Patient Category AD

FMP 20 SSN 111-22-3333

Facility/Unit USNH Naples Branch Navt Rank/Grade 0-3

Civilian Occupation French Fry Cook Religion No Preference Race Caucasian

Home Phone Work Phone DSN 629-6382

Address USNH Naples

PSC 827

City FPO State AE Zip 09617

Sponsor's name Relationship

Complete PT Demographics then click on "Ophthalmology"

Study Data

Accession Number 1127401615 Modality SC

Study description EYE-SL

Operator Mick Jager

Physician Data




Referring Physician Dr. Evil Referring Physician Phone DSN 629-6382

Enter the appropriate study description:

- EYE-IVFA
- EYE-FUNDUS
- EYE-HVF
- EYE-OCT
- EYE SL

You can search by modality, so note the modality type is "SC"

Import from Badge file... Import from Insurance Card file... Print label...

  ☐ Publish into CD ☐ Print  Hint: Use (CTRL+) right mouse button to set image series, SHIFT+ right mouse button to rotate, ...

Immunographics | Dermatology | PATH | Ophthalmology | Images

Patient Medical Information | Patient Medical Information | Segment Evaluation | Provisional Diagnosis

Allergies

Medications
Rx
OTC
Herb/Diet/Vitamin Supp

Chief Complaint

Pain
☒ Yes ☐ No Location
Character Intensity
Oriented to time/place and person ☒ Yes ☐ No Mood/affect normal ☒ Yes ☐ No

Dist VA
cc/sc

Near VA
cc/sc

Current Rx
OD ☐ NCT OD
OS ☒ Goldmann OS
☐ Other Time

Dilation
Time ☐ 0.5% proparacaine ☒ 2.5% phenylephrine
☒ 1.0% tropicamide ☐ 1.0% cyclopentolate

Fill out exam
information by
clicking on tabs

☐ Publish into CD☐ Print

Settings...

Current series (1-100)

1



Hint: Use (CTRL+) right mouse button
to set image series, SHIFT+ right
mouse button to rotate, ...

Ophthalmographics | Dermatology | PATH | Ophthalmology | Images

Patient Medical Information | Patient Medical Information | Segment Evaluation | Provisional Diagnosis

Cover Test

D Ortho N Ortho

Pupils PERRL (-) APD

Conf VF Full to finger count

Motility FULL no restrictions

Subjective

OD PH 20/ NI

OS 20/

Final Rx

OD

OS Add

Additional Tests

BP 120 / 80 Pulse 67 Temp 97.5

Color 9/9 PIP

☐ Publish into CD☐ Print [Settings...](#)

Current series (1-100)

1



Hint: Use (CTRL+) right mouse button
to set image series, SHIFT+ right
mouse button to rotate, ...

[Ophthalmographics](#) | [Dermatology](#) | [PATH](#) | [Ophthalmology](#) | [Images](#)[Patient Medical Information](#) | [Patient Medical Information](#) | [Segment Evaluation](#) | [Provisional Diagnosis](#)

Anterior Segment Evaluation

Lids/Lashes	NL
Conjunctiva	2+ Inject
Cornea	Clear
Anterior Chamber VH: 1 2 3 4	3
Iris	CL
Lens	2 + Cell and Flare

Posterior Segment Evaluation

☐ Undilated ☒ Dilated
☐ Direct ☒ BIO ☐ 90D ☐ 78D ☐ 66D ☐ 60D

OD		C/D	.3/.3 OU
OS			
OU			

Vessels	NL
Macula	NL
Vitreous	NL
Periphery	Atrophic hole



Publish into CD



Print

Settings...

Current series (1-100)

1



Hint: Use (CTRL+) right mouse button
to set image series, SHIFT+ right
mouse button to rotate, ...

Immunographics

Dermatology

PATH

Ophthalmology

Images

Patient Medical Information

Patient Medical Information

Segment Evaluation

Provisional Diagnosis

Provisional Diagnosis, Questions, Comments

Please enter provisional diagnosis in the space below

Badness on or about the eye

Specify questions from provider or patient to Ophthalmologist

Medivac?

Additional relevant information

When you are finished
entering exam data
click on "Images"



☐ Publish into CD

☐ Print

Settings...

Current series (1-100)

1



Hint: Use (CTRL+) right mouse button to set image series, SHIFT+ right mouse button to rotate, ...

Imagographics

Dermatology

PATH

Ophthalmology

Images

Drop images HERE from capture devices in the lower pane...

Scanner (TWAIN)

VIDAR

DENTAL DR

ECG

Spiro

FILE...

Select Source

Acquire

Create from file...

TWAIN Data Source Manager Open

Click on source

☐ Publish into CD☐ Print

Settings...

Current series (1-100)

1



Hint: Use (CTRL+) right mouse button to set image series, SHIFT+ right mouse button to rotate, ...

Imaging | Dermatology | PATH | Ophthalmology | Images


Drop images HERE from capture devices in the lower pane...

Scanner (TWAIN) | VIDAR | DENTAL DR

Select Source

Acqu

TWAIN Data Source Manager Open

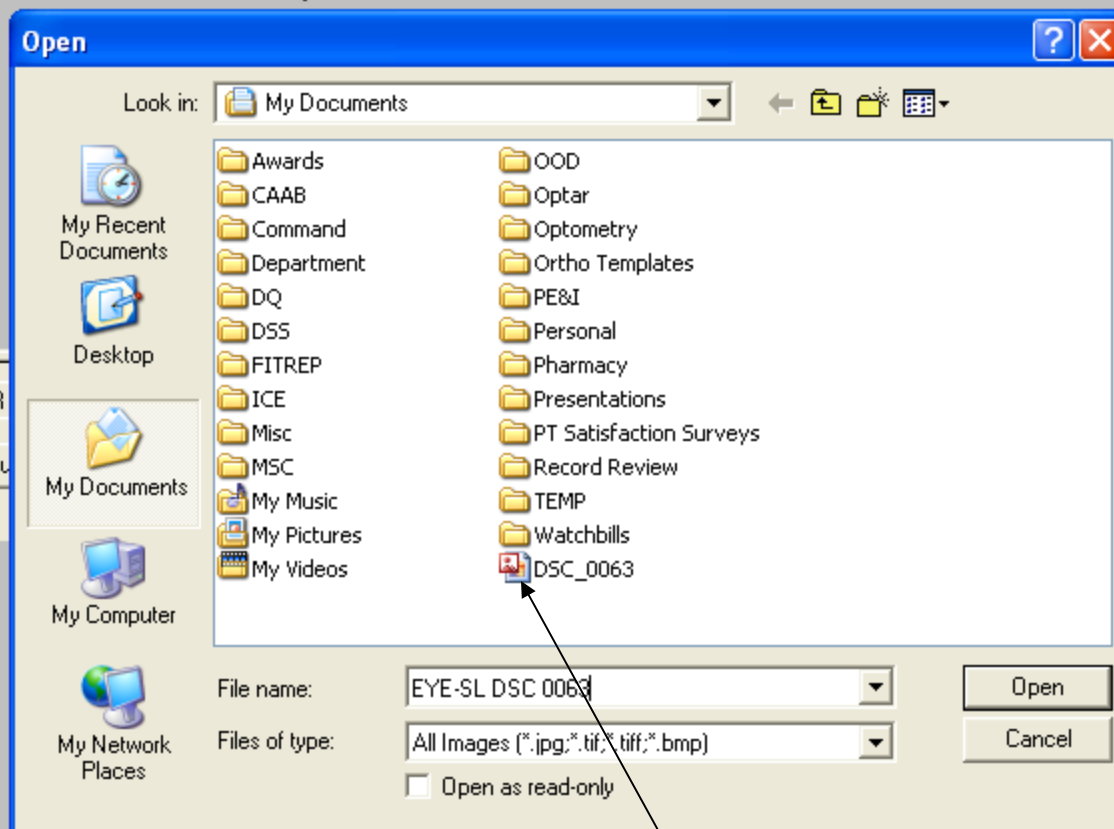

Select
File



☐ Publish into CD

☐ Print

Settings...

Current series (1-100)

1



Hint: Use (CTRL+) right mouse button to set image series, SHIFT+ right mouse button to rotate, ...

Imagographics

Dermatology

PATH

Ophthalmology

Images

Drop images HERE from capture devices in the lower pane...

Scanner (TWAIN)

VIDAR

DENTAL DR

ECG

Spiro

FILE...


Select Source

Acquire

Create from file...

TWAIN Data Source Manager Open

: NONE (COLOR)

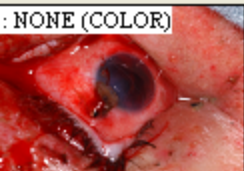

Drag file to top window

☐ Publish into CD☐ Print

Settings...

Current series (1-100)

1



Hint: Use (CTRL+) right mouse button to set image series, SHIFT+ right mouse button to rotate, ...

omographics

Dermatology

PATH

Ophthalmology

Images

1 (COLOR)



Click on "upload Image"
Icon

canner (TWAIN)

VIDAR

DENTAL DR

ECG

Spiro

FILE...

Select Source

Acquire

Create from file...

TWAIN Data Source Manager Open

: NONE (COLOR)





☐ Publish into CD

☐ Print

Settings...

Current series (1-100)

1



Hint: Use (CTRL+) right mouse button to set image series, SHIFT+ right mouse button to rotate, ...

Imagographics

Dermatology

PATH

Ophthalmology

Images

1 (COLOR)



Scanner (TWAIN) | VIDAR | DENTAL DR | ECG | Spine

Select Source

Acquire

TWAIN Data Source Manager Open

: NONE (COLOR)



Medweb Capture Application



Are you sure you want to upload all selected images and new reports?

OK

Cancel

Click "OK"

☐ Publish into CD☐ Print [Settings...](#)

Current series (1-100)


1



Hint: Use (CTRL+) right mouse button to set image series, SHIFT+ right mouse button to rotate, ...

Imaging | Dermatology | PATH | Ophthalmology | Images

1 (COLOR)

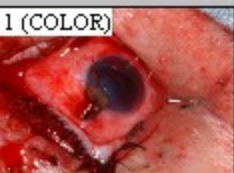

Scanner (TWAIN) | VIDAR | DENTAL DR | ECG | Spiro

Select Source

Acquire

TWAIN Data Source Manager Open

: NONE (COLOR)

**Security Alert**

Information you exchange with this site cannot be viewed or changed by others. However, there is a problem with the site's security certificate.



The security certificate was issued by a company you have not chosen to trust. View the certificate to determine whether you want to trust the certifying authority.



The security certificate has expired or is not yet valid.



The name on the security certificate is invalid or does not match the name of the site

Do you want to proceed?

Yes

No

View Certificate

Click
"Yes"



☐ Publish into CD

☐ Print

Settings...

Current series (1-100)

1



Hint: Use (CTRL+) right mouse button to set image series, SHIFT+ right mouse button to rotate, ...

Imaging | Dermatology | PATH | Ophthalmology | Images

1 (COLOR)



Scanner (TWAIN) | VIDAR | DENTAL DR | ECG | Spiro | FILE

Select Source

Acquire

Create

TWAIN Data Source Manager Open

: NONE (COLOR)



Connect to 131.158.69.236



Patient Record Access

User name:

pgunther

Password:

.....

OK

Cancel

Enter user name and password

☐ Publish into CD☐ Print [Settings...](#)

Current series (1-100)


1



Hint: Use (CTRL+) right mouse button to set image series, SHIFT+ right mouse button to rotate, ...

Imagographics | Dermatology | PATH | Ophthalmology | Images

1 (COLOR)

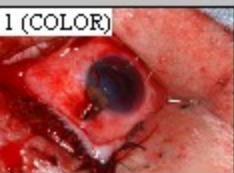

Scanner (TWAIN) | VIDAR | DENTAL DR | ECG | Spiro

Select Source

Acquire

TWAIN Data Source Manager Open

: NONE (COLOR)

**Security Alert**

Information you exchange with this site cannot be viewed or changed by others. However, there is a problem with the site's security certificate.



The security certificate was issued by a company you have not chosen to trust. View the certificate to determine whether you want to trust the certifying authority.



The security certificate has expired or is not yet valid.



The name on the security certificate is invalid or does not match the name of the site

Do you want to proceed?

Yes

No

View Certificate

Click "View Certificate"

☐ Publish into CD☐ Print Settings...

Current series (1-100)


1



Hint: Use (CTRL+) right mouse button
to set image series, SHIFT+ right
mouse button to rotate, ...

Imaging | Dermatology | PATH | Ophthalmology | Images

1 (COLOR)

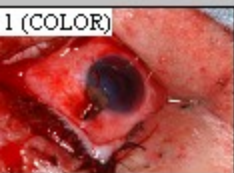

Scanner (TWAIN) | VIDAR | DENTAL DR | ECG | Spiro

Select Source

Acquire

TWAIN Data Source Manager Open

: NONE (COLOR)



Click "OK" to install
Certificate

Certificate

General Details Certification Path



Certificate Information

**This certificate cannot be verified up to a trusted
certification authority.**

Issued to: devel**Issued by:** development**Valid from** 12/14/2001 **to** 9/9/2004

Install Certificate...

Issuer Statement

OK

☐ Publish into CD☐ Print

Settings...

Current series (1-100)

1



Hint: Use (CTRL+) right mouse button
to set image series, SHIFT+ right
mouse button to rotate, ...

Imaging


Dermatology

PATH

Ophthalmology

Images

1 (COLOR)

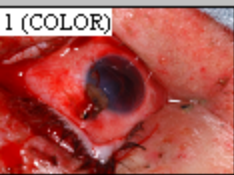

Scanner (TWAIN)

VIDAR

DENTAL DR

ECG


Spiro

Select Source

Acquire

TWAIN Data Source Manager Open

: NONE (COLOR)

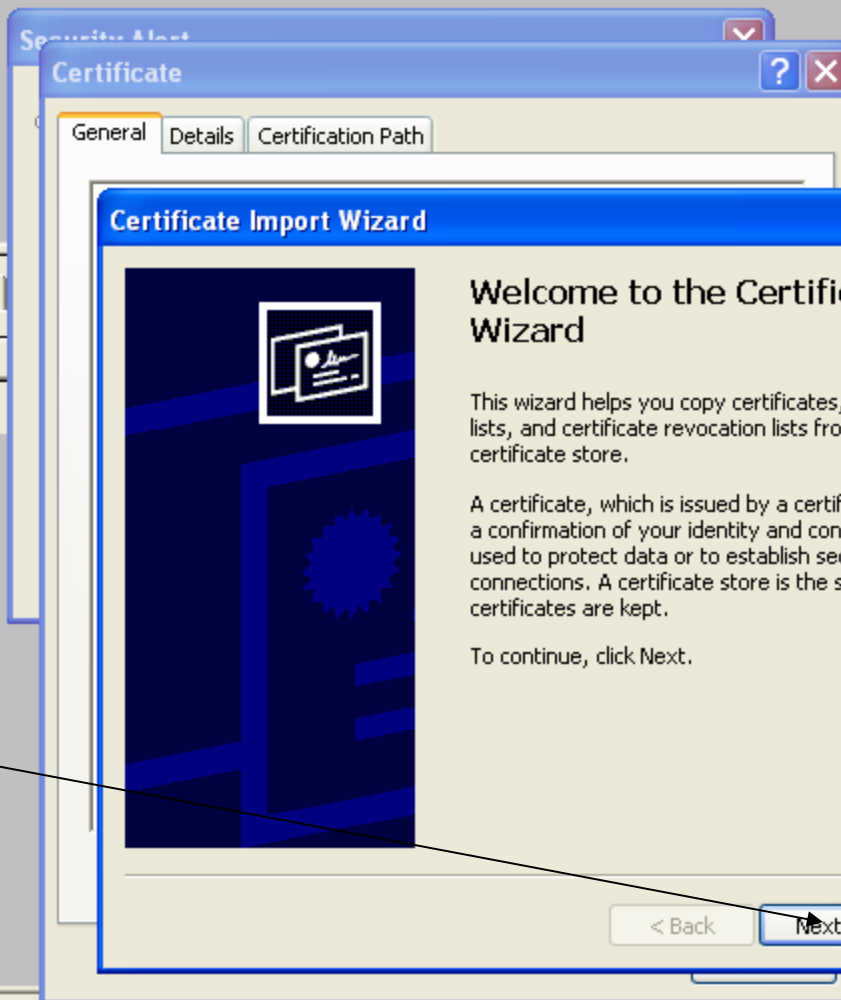

Certificate Import Wizard

Welcome to the Certificate Import Wizard

This wizard helps you copy certificates, certificate trust lists, and certificate revocation lists from your disk to a certificate store.

A certificate, which is issued by a certification authority, is a confirmation of your identity and contains information used to protect data or to establish secure network connections. A certificate store is the system area where certificates are kept.

To continue, click Next.

< Back

Next >

Cancel

Complete Import
Wizard

☐ Publish into CD☐ Print

Settings...

Current series (1-100)

1



Hint: Use (CTRL+) right mouse button
to set image series, SHIFT+ right
mouse button to rotate, ...

Imaging | Dermatology | PATH | Ophthalmology | Images

1 (COLOR)



Scanner (TWAIN) | VIDAR | DENTAL DR | ECG | Spiro

Select Source

Acquire

TWAIN Data Source Manager Open

: NONE (COLOR)



Certificate

General

Details

Certification Path

Certificate Import Wizard

Certificate Store

Certificate stores are system areas where certificates are kept.

Windows can automatically select a certificate store, or you can specify a location for

- ☒ Automatically select the certificate store based on the type of certificate
- ☐ Place all certificates in the following store

Certificate store:

Browse...

< Back

Next >

Cancel

☐ Publish into CD☐ Print

Settings...

Current series (1-100)

1



Hint: Use (CTRL+) right mouse button
to set image series, SHIFT+ right
mouse button to rotate, ...

mographics | Dermatology | PATH | Ophthalmology | Images

1 (COLOR)



canner (TWAIN) | VIDAR | DENTAL DR | ECG | Spiro

Select Source

Acquire

TWAIN Data Source Manager Open

: NONE (COLOR)



Security Alert


Certificate

General

Details

Certification Path

Certificate Import Wizard

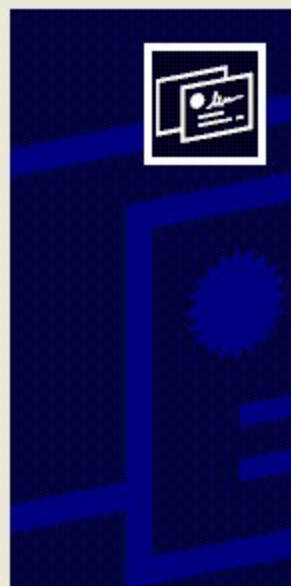

Completing the Certificate Import Wizard

You have successfully completed the Certificate Import wizard.

You have specified the following settings:

Certificate Store Selected	Automatically determined by t
Content	Certificate

< Back

Finish

Cancel

☐ Publish into CD☐ Print

Settings...

Current series (1-100)


1



Hint: Use (CTRL+) right mouse button to set image series, SHIFT+ right mouse button to rotate, ...

Immunographics | Dermatology | PATH | Ophthalmology | Images

1 (COLOR)

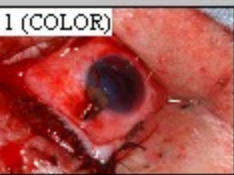


Scanner (TWAIN) | VIDAR | DENTAL DR | ECG | Spiro

Select Source

Acquire

TWAIN Data Source Manager Open

: NONE (COLOR)

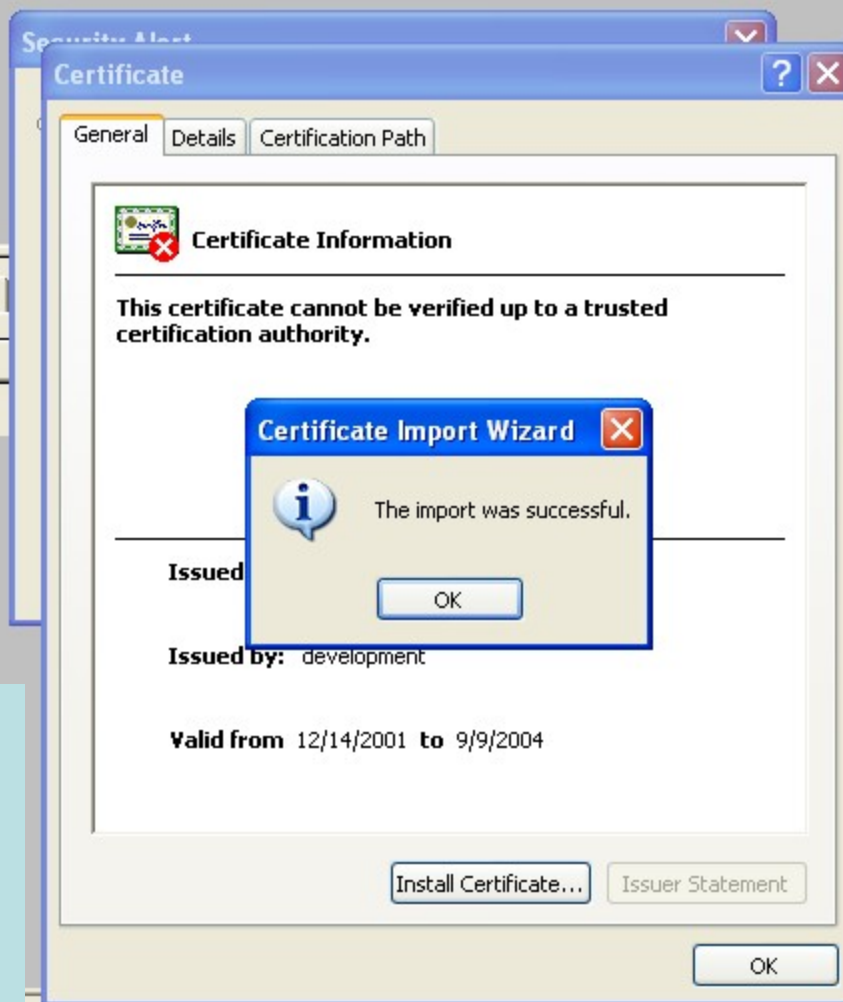

The image will upload automatically and you are done! Now email consulting provider so they know to check the Medweb.

- Home
- New Patients
- Patient schedule
- Machine schedule
- Monitor
- Personal Profile
- Statistics
- Logout
- Administrator Site
- Admin Site SSL

NNMC Bethesda

The consultant will email the referring provider letting them know that the consult is complete and posted on Medweb. The referring provider can then log back into Bethesda's server to view the results.

Announcements: 501.293.3812



Local site info

Address and driving directions

Address: United States
8901 Wisconsin Ave, Bethesda MD

- medweb
- Home
- New Patients
- Patient schedule
- Machine schedule
- Monitor
- Personal Profile
- Statistics
- Logout
- Administrator Site
- Admin Site SSL

NNMC Bethesda

[Ophthalmologist] Medweb tele-medicine intranet software - Version 4.5.20 - Last updated 2005-08-0

New Study On Call Preferences Save settings View Changes Manuals

Get latest plugin

Your license has 3376 days remaining.
Your Medweb Server manager can update the license.

Limit:	Status:	Find:	Modality:	In last:	Group By:	Query:
10	Unread		SC	(any)	Site/Station	Local

Search:	Advanced	Current	Clear	User: pgunther	Drive status: 110033M available for studies
				0 images left to send	0 images left to compress

Result Pages: 1 2 [Next>]

☐ View only pre-selected images

	Name	ID	Type	Imgs	Date	Time	Mgmt
Site: Nnmc Bethesda Station name:NNMC Bethesda Modality:SC Compression: Low							
<input type="checkbox"/>	Gunther Peter	1127401615	EYE-SL	1	2005-09-22	10:48:13	
<input type="checkbox"/>	Lynch James	082660084	eye-fundus	1	2005-09-22	06:21:25	
<input type="checkbox"/>	Case Test	1127227696	Test Case	1	2005-09-20	11:33:47	
<input type="checkbox"/>	Blice Jeff			4	2005-09-19	12:42:15	
Site: Okinawa Compression: Low							
<input type="checkbox"/>	Test Hancock			1	2005-07-15	03:27:41	
Site: Rota Station name:Rota Medweb Modality:SC Compression: Low							
<input type="checkbox"/>	Lynch James	082660084	EYE-TOPO	2	2005-09-21	08:19:03	
<input type="checkbox"/>	Rota Optometry Testimage	1127227696		1	2005-09-20	16:49:44	

Click on the icon to view the completed consult

Click here to view or modify Web-reports for this study

Pt.: Gunther Peter

DESCRIPTION EYE-SL

History and
Provisional diagnosis**Study reports:****Add**

Created	Description	Format	Mgmt
2005-09-23 11:45:59	2005/09/23.11:45:59 - Jeffrey P. Blice: EYE-SL SIGNED BY: Jeffrey P. Blice, NMMC Bethesda, 301-295-0373, JPBlice (2005/09/23.11:46:00) Report status: FINAL RESULTS; ORDER IS COMPLETE AND VERIFIED	html	
2005-09-22 11:00:00	2005/09/22 11:00:00 - : EYE-SL SIGNED BY: , NMMC Bethesda (2005/09/22 11:00:00) Report status: FINAL RESULTS; ORDER IS COMPLETE AND VERIFIED	html	
2005-09-22 10:48:09	OPHTHALMOLOGY: ProvDiag	text/html	
2005-09-22 10:48:03	OPHTHALMOLOGY: SegEval	text/html	
2005-09-22 10:47:20	OPHTHALMOLOGY: MedInfo2	text/html	
2005-09-22 10:47:13	OPHTHALMOLOGY: MedInfo1	text/html	
2005-09-22 10:44:44	DEMOGRAPHICS: Patient Information	application/x-medweb-alist-multiline	

**Click on the signed
consult**

Mgmt

Compression: Low



Compression: Low



Compression: Low



